



"Tomorrow's Instruments... Today"™

TANNAS *QUANTUM* OXIDATION TESTER

Oxidation Stability – RPVOT & TFOUT

- ASTM D 2272, IP 229, D 4742, D 7098 – used with new & in-service oils, steam turbine oils, insulating oils and gasoline engine oils.
- Only non-Liquid RPVOT & TFOUT oxidation bath available
- New direct 'Dry-Cylinder' sample heating -- eliminates hot, hazardous, liquid bath mess and odor
- Rapid turn-around in test capabilities due to independent sample testing -- estimated to double productivity with multi-unit setup
- Very small, bench-top footprint
- Does not require placement in hood -- simple venting of odors through plastic tubing to scrubber or vent
- Quad-Automation Package with laptop, monitors & records up to four *Quantum* units at the same time or separately



Hot oil baths of varying sizes have long been used in the industry for running the Rotating Pressure Vessel Oxidation Test or RPVOT (formerly known as RBOT) and the Thin-Film Oxygen Uptake test (TFOUT). These bath types expose the operator to hot oil fumes, are prone to leak, need regular maintenance and are generally messy and potentially hazardous to operate.

As laboratory facilities are under increasing scrutiny to improve working conditions and encourage healthy work environments, alternative approaches to running these tests is desirable.

In 2003, Tannas Co. set out to effectively eliminate the hazards and mess of hot oil baths and has recently completed the development of the Tannas *Quantum* Oxidation Tester, **the only non-liquid 'dry-cylinder' heating system available for running these critical tests.**

TANNAS CO.



4800 James Savage Rd.
Midland, MI 48642 USA
Phone: 989-496-2309
Fax: 989-496-3438
Website: www.savantgroup.com
Email: tannas@savantgroup.com

Tannas *Quantum* Oxidation Tester

Dimensions	Bench-top 8"(w) x 15"(d) x 13" (20 x 38 x 33 cm)
Weight	~20 lbs. (9 kg)
Voltage	120/220 VAC, Single Phase, 15 amp. 50/60 Hz
Heating Medium	'Dry-Cylinder' heating system – no hot oil bath
Testing Capacity	Single position stainless steel pressure vessel w/pressure transducer Designed for multi-unit alignment, each independent
Test Parameter Capabilities	Temperature – 200°C ±0.1°C recommended max. Oxygen Charge – 100 ±0.1 psi recommended max. Vessel Rotation – Variable speed control
Output	Continuous temperature & oxygen pressure readout -Digital USB to Laptop -Dsub9 to Analog strip-chart recording
Safety	Current limiting fuses Over-pressure sensor & relief Over-temperature Cut-off Fuse



Data taken from four *Quantum* instrument qualifications show close comparison to that of the traditional liquid bath on a low and relatively high reference oil. *Quantum* results also show good precision.

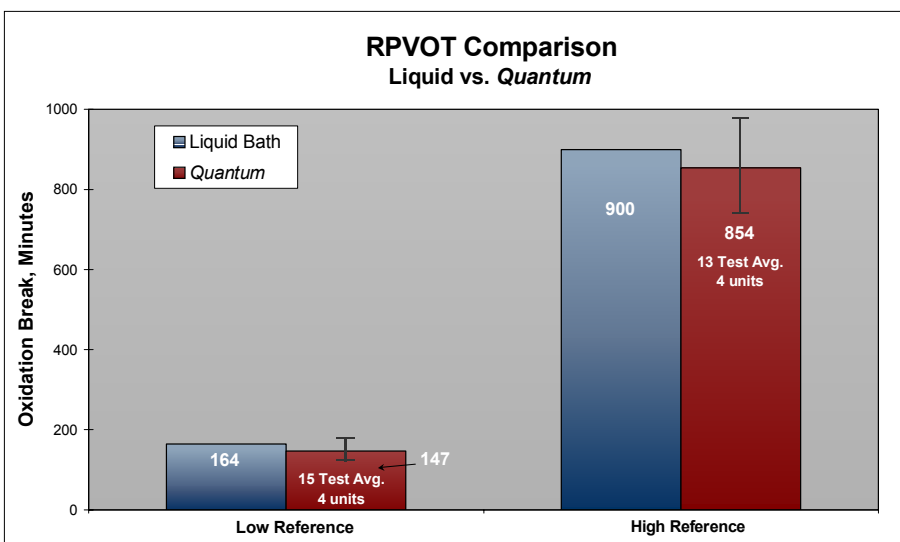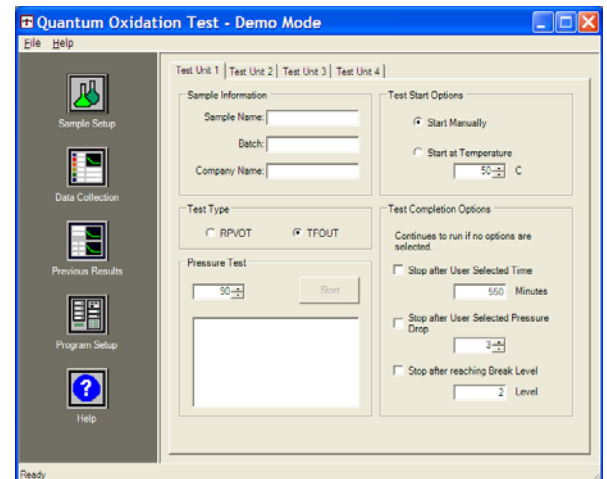

The *Quantum* software package provides automated pressure and temperature monitoring and end point detection for up to four independent test units. End of test criteria can be selected by Time Duration, Pressure Drop or Break level. DAQ acquisition with USB interface to included laptop...
Many more features, call for details.



TANNAS CO. ◊ 4800 James Savage Rd. ◊ Midland, MI ◊ 48642 ◊ USA
989-496-2309 ◊ Fax: 989-496-3438 ◊ Website: www.savantgroup.com/tannas.sht
Email: tannas@savantgroup.com